



UROVISION • UROTECH

A MEDI-GLOBE COMPANY

Pediatric Urology Product Catalogue



Your Partner in
Endourology

CONTENTS

WE CARE FOR PEDIATRIC UROLOGY

There are many reasons why urological examinations and interventions are necessary in children and adolescents. Recurring urinary tract infections, malformations in the urogenital tract, hypospadias or oncological diseases of the urogenital tract are only a few indications that make interventions necessary. As a whole, urological diseases in childhood are not uncommon.

As a leading provider of single-use medical products for endourological interventions, WE for sure provide a special range of products tailored to the needs of Pediatric Urology. Besides standard products in smaller sizes we have specially developed, innovative endourological products to support therapy.

UROTECH/UROVISION - Your Partner in Endourology.
 With innovative solutions that facilitate interventions and spare children discomfort and pain as much as possible.

WE ARE URO!



STONE BASKETS TIPLESS
[Page 4](#)



URETERAL CATHETERS
[Page 5](#)



URETERAL STENTS
[Page 6](#)



NEPHROSTOMY
[Page 16](#)




SUPRAPUBIC BLADDER DRAINAGE
[Page 20](#)



TRANSURETHRAL BLADDER DRAINAGE
[Page 24](#)



SPECIAL PRODUCTS
[Page 26](#)

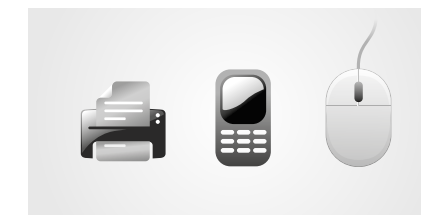

ORDERING INFORMATION
[Back Page](#)



STONE BASKETS TIPLESS


STONE BASKETS
High elastic Nitinol material

Handle acc. to Lahme for all stone baskets

- De- and reconnectable handle:
Endoscope can be removed while basket is situated. Handle and basket can get reattached afterwards
- Fixation of the stone in the basket due to spring mechanism
- Ergonomic handle design with non-slip surface
- Basket opens and closes within image area:
No adjustment of endoscope necessary



„TIPLESS“

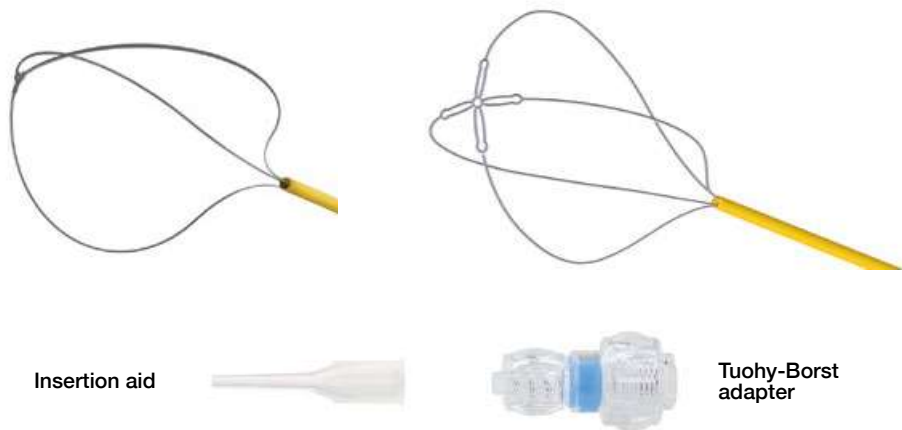


- "Smallest basket in the world" with 1.2 FR/CH
- Best rinse flow characteristics and visibility
- For picking up the smallest stone concretions even at the edge of the tissue
- With Tuohy-Borst adapter for fixation to the working channel of the endoscope
- Insertion aid protects the tip of the basket during insertion into the endoscope

Article code	Size (FR/CH)	Wires	Working Length (cm)	Basket Length (mm)	Basket Ø (mm)	Basket-form
STN-T123-120	1.2	3	120	12	8	straight
STN-T184-120	1.8	4	120	15	11	straight

Packaging unit = 3 Pcs

URETERAL CATHETER

URETERAL CATHETER
For retrograde Urography

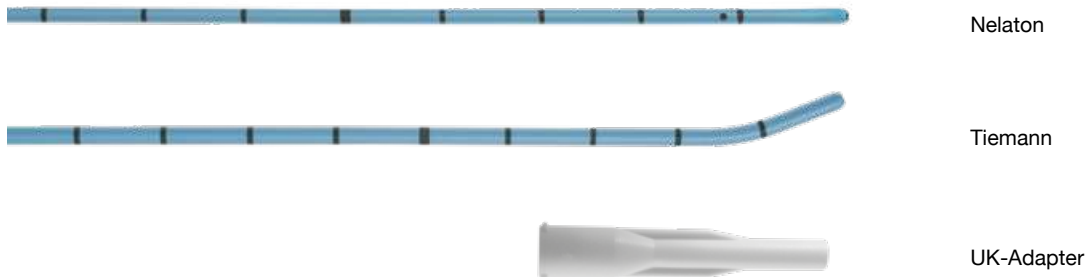

- With universal adapter for injection of contrast medium
- Good X-ray visibility
- With inward stainless-steel stylet

Nelaton-tip				Packaging unit = 10 Pcs
Article code closed tip	Article code open tip	Size (CH/FR)	Length (cm)	Max Guide wire (inch)
UK-500370	UK-600370	3	70	0.025"
UK-500470	UK-600470	4	70	0.035"


Tiemann-tip				Packaging unit = 10 Pcs
Article code closed tip	Article code open tip	Size (CH/FR)	Length (cm)	Max Guide wire (inch)
UK-510370	UK-540370	3	70	0.025"
UK-510470	UK-540470	4	70	0.035"

URETERAL STENTS

YELLOW-STAR Double-J Ureteral Stent-Set
PC-treated



SET COMPONENTS

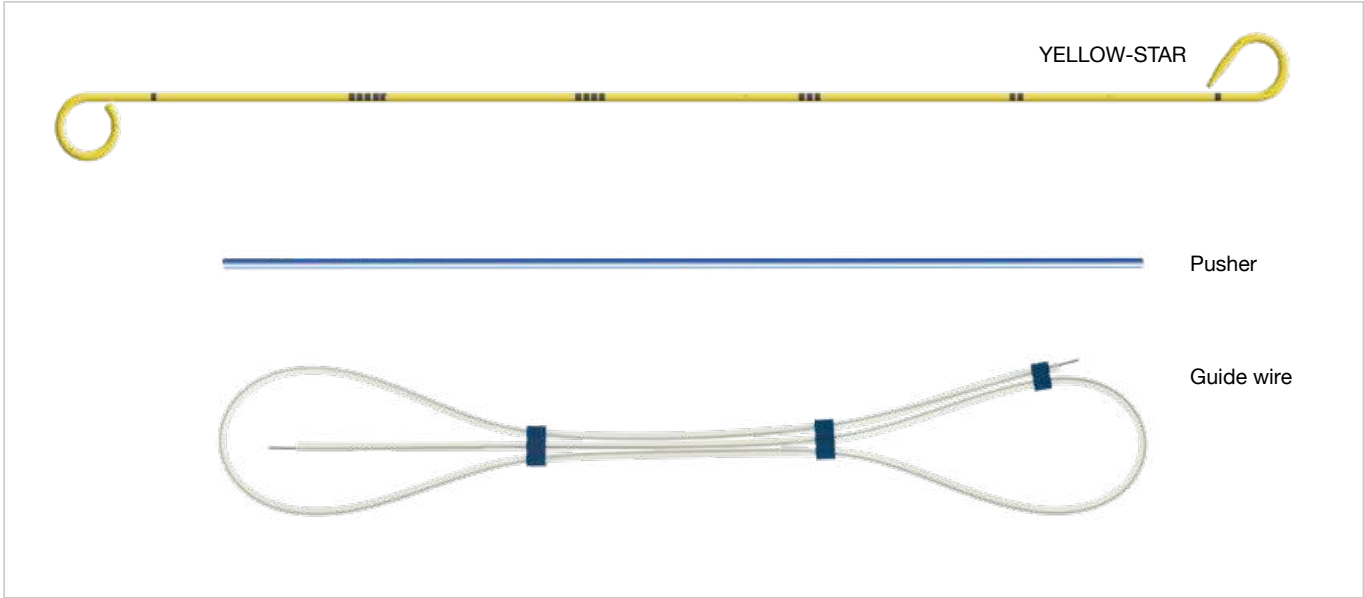

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ PC-treatment to minimize risk of infection and encrustation, resulting in extended indwelling time
- ◆ Drainage openings in pigtails and shaft
- ◆ Open tip
- ◆ Markings for exact positioning
- ◆ X-ray visible

Pusher

- ◆ Length: 40 cm

Guide wire standard

- ◆ PTFE-coated

URETERAL STENTS

Open tip							Packaging unit = 10 Pcs
Article code w/o Guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
EP-414710	EP-364710	4.7	10	125	0.028"	straight	standard
EP-414712	EP-364712	4.7	12	125	0.028"	straight	standard
EP-414715	EP-364715	4.7	15	125	0.028"	straight	standard
EP-414718	EP-364718	4.7	18	125	0.028"	straight	standard
EP-414720	EP-364720	4.7	20	125	0.028"	straight	standard
EP-414722	EP-364722	4.7	22	125	0.028"	straight	standard

Guide wire standard – PTFE-coated stainless steel						Packaging unit = 10 Pcs
Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft	
ST-291528	0.028"	150	straight	standard	standard	

UROLIX guide wire "hydrophilic"						Packaging unit = 10 Pcs
Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft	
UXN-SS-2515	0.025"	150	straight	standard	standard	

Platin Pusher				Packaging unit = 10 Pcs
Article code	Size (FR/CH)	Length (cm)	Max Guide wire Ø	
PU-0545-RM	5	45	0.028"	
PU-0570-RM	5	70	0.028"	


URETERAL STENTS

PREMIUM YELLOW STAR PLATIN PUSHER-LINE

For visually secured placement of the ureteral stent



SET COMPONENTS

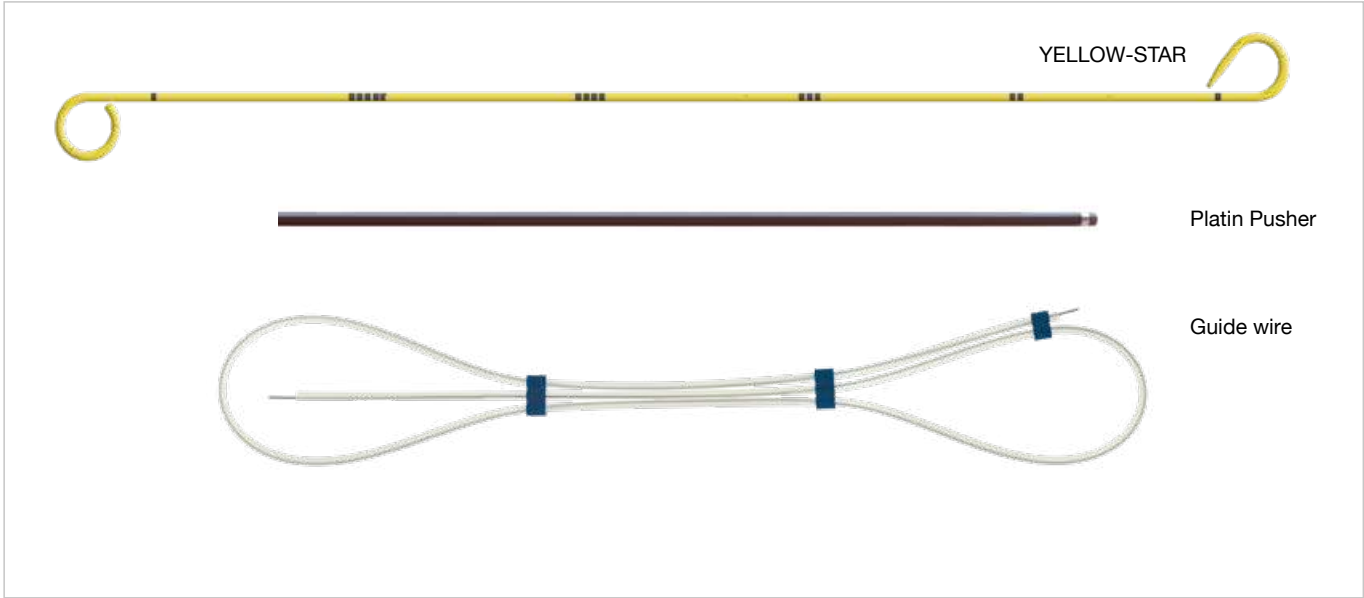

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ PC-treatment to minimize risk of infection and encrustation, resulting in extended indwelling time
- ◆ Drainage openings in pigtails and shaft
- ◆ Open tip
- ◆ Markings for exact positioning
- ◆ X-ray visible

Platin Pusher

- ◆ Improved radiopaque visibility due to platinum ring at the tip
- ◆ Increased safety during retrograde insertion into the bladder

Guide wire standard

- ◆ PTFE-coated

URETERAL STENTS

Open tip							Packaging unit = 10 Pcs
Article code w/o Guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
-	YSPL-364708	4.7	8	125	0.028"	straight	standard
YSPL-414710	YSPL-364710	4.7	10	125	0.028"	straight	standard
-	YSPL-364710-F	4.7	10	125	0.028"	straight	standard
YSPL-414712	YSPL-364712	4.7	12	125	0.028"	straight	standard
-	YSPL-364712-F	4.7	12	125	0.028"	straight	standard
YSPL-414715	YSPL-364715	4.7	15	125	0.028"	straight	standard
YSPL-414715-F	YSPL-364715-F	4.7	15	125	0.028"	straight	standard
YSPL-414718	YSPL-364718	4.7	18	125	0.028"	straight	standard
YSPL-414718-F	YSPL-364718-F	4.7	18	125	0.028"	straight	standard
YSPL-414720	YSPL-364720	4.7	20	125	0.028"	straight	standard
YSPL-414720-F	YSPL-364720-F	4.7	20	125	0.028"	straight	standard
YSPL-414722	YSPL-364722	4.7	22	125	0.028"	straight	standard
YSPL-414722-F	YSPL-364722-F	4.7	22	125	0.028"	straight	standard

Guide wire standard - PTFE-coated stainless steel						Packaging unit = 10 Pcs
Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft	
ST-291528	0.028"	150	straight	standard	standard	

UROLIX Guide wire "hydrophilic"						Packaging unit = 10 Pcs
Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft	
UXN-SS-2515	0.025"	150	straight	standard	standard	


Platin Pusher				Packaging unit = 10 Pcs
Article code	Size (FR/CH)	Length (cm)	Max Guide wire Ø	
PU-0545-RM	5	45	0.028"	
PU-0570-RM	5	70	0.028"	

Variants “-F” are with string

URETERAL STENTS

WHITE-STAR Double-J Ureteral Stent-Set VISIOSTAR steerable set variant

SET COMPONENTS

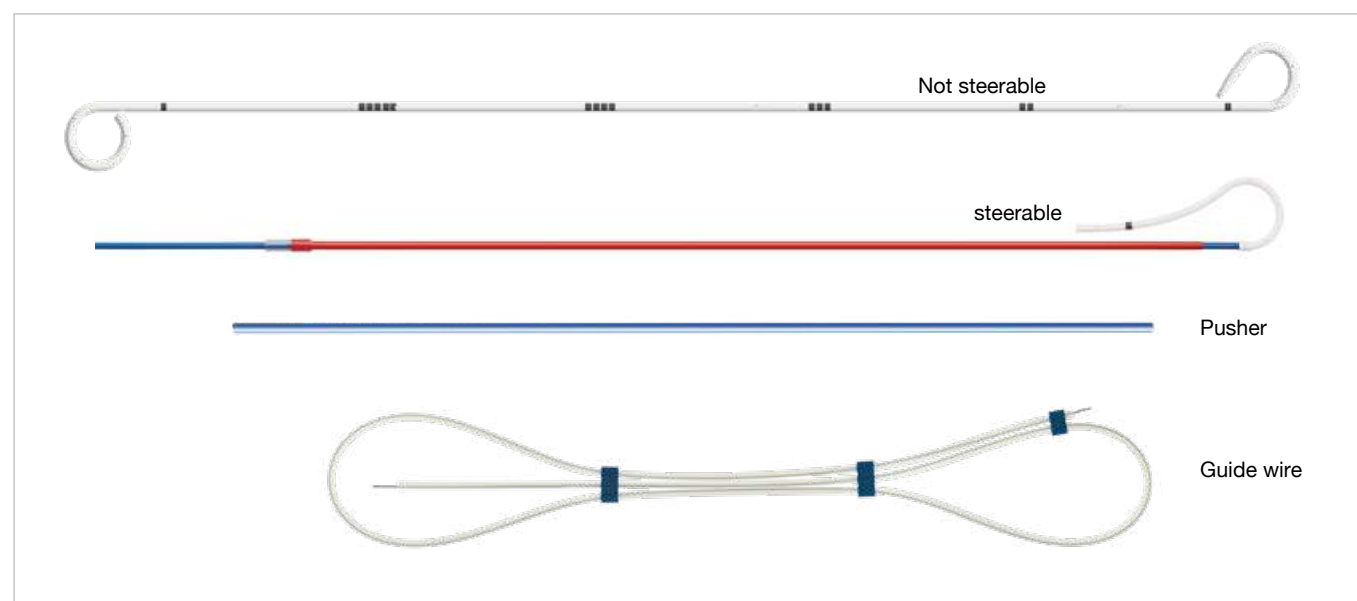

Ureteral Stent

- ♦ Polyurethane (PUR)
- ♦ Drainage openings in pigtails and shaft
- ♦ Open tip
- ♦ Markings for exact positioning
- ♦ X-ray visible

Pusher

- ♦ Length: 40 cm

Guide wire standard

- ♦ PTFE-coated

VISIOSTAR steerable

- ♦ Catheter can be easily rotated and moved up or down
- ♦ Secure positioning of the ureteral stent
- ♦ Pre-assembled complete set

URETERAL STENTS

Open tip not steerable

Packaging unit = 10 Pcs

Article code w/o guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
ST-197308	ST-190308	3	8	150	0.018"	straight	standard
ST-197310	ST-190310	3	10	150	0.018"	straight	standard
ST-197312	ST-190312	3	12	150	0.018"	straight	standard
ST-197315	ST-190315	3	15	150	0.018"	straight	standard
ST-197318	ST-190318	3	18	150	0.018"	straight	standard
ST-197320	ST-190320	3	20	150	0.018"	straight	standard
-	ST-190322	3	22	150	0.018"	straight	standard
EN-410408	EN-360408	4	08	125	0.025"	straight	standard
EN-410410	EN-360410	4	10	125	0.025"	straight	standard
EN-410412	EN-360412	4	12	125	0.025"	straight	standard
EN-410415	EN-360415	4	15	125	0.025"	straight	standard
EN-410418	EN-360418	4	18	125	0.025"	straight	standard
EN-410420	EN-360420	4	20	125	0.025"	straight	standard
-	EN-360422	4	22	125	0.025"	straight	standard
EN-410508	-	4.7	8	125	0.028"	straight	standard
EN-410510	EN-360510	4.7	10	125	0.028"	straight	standard
EN-410512	EN-360512	4.7	12	125	0.028"	straight	standard
EN-410515	EN-360515	4.7	15	125	0.028"	straight	standard
EN-410518	EN-360518	4.7	18	125	0.028"	straight	standard
EN-410520	EN-360520	4.7	20	125	0.028"	straight	standard
EN-410522	EN-360522	4.7	22	125	0.028"	straight	standard

Open tip steerable

Packaging unit = 10 Pcs

Article code with guide wire	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
US-110508	4.7	8	150	0.028"	straight	standard
US-110510	4.7	10	150	0.028"	straight	standard
US-110512	4.7	12	150	0.028"	straight	standard
US-110515	4.7	15	150	0.028"	straight	standard
US-110518	4.7	18	150	0.028"	straight	standard
US-110520	4.7	20	150	0.028"	straight	standard
US-110522	4.7	22	150	0.028"	straight	standard

Guide wire standard - PTFE-coated stainless steel

Packaging unit = 10 Pcs

Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft
ST-291518	0.018"	150	straight	standard	standard
ST-291528	0.028"	150	straight	standard	standard

UROLIX guide wire "hydrophilic"


Packaging unit = 10 Pcs

Article code	Ø (inch)	Length (cm)	Tip	Stiffness Tip	Stiffness Shaft
UXN-SS-2515	0.025"	150	straight	standard	standard

URETERAL STENTS

MAGNETIC BLACK-STAR® Double-J Ureteral Stent-Set and Retrieval Device with magnet

SET COMPONENTS

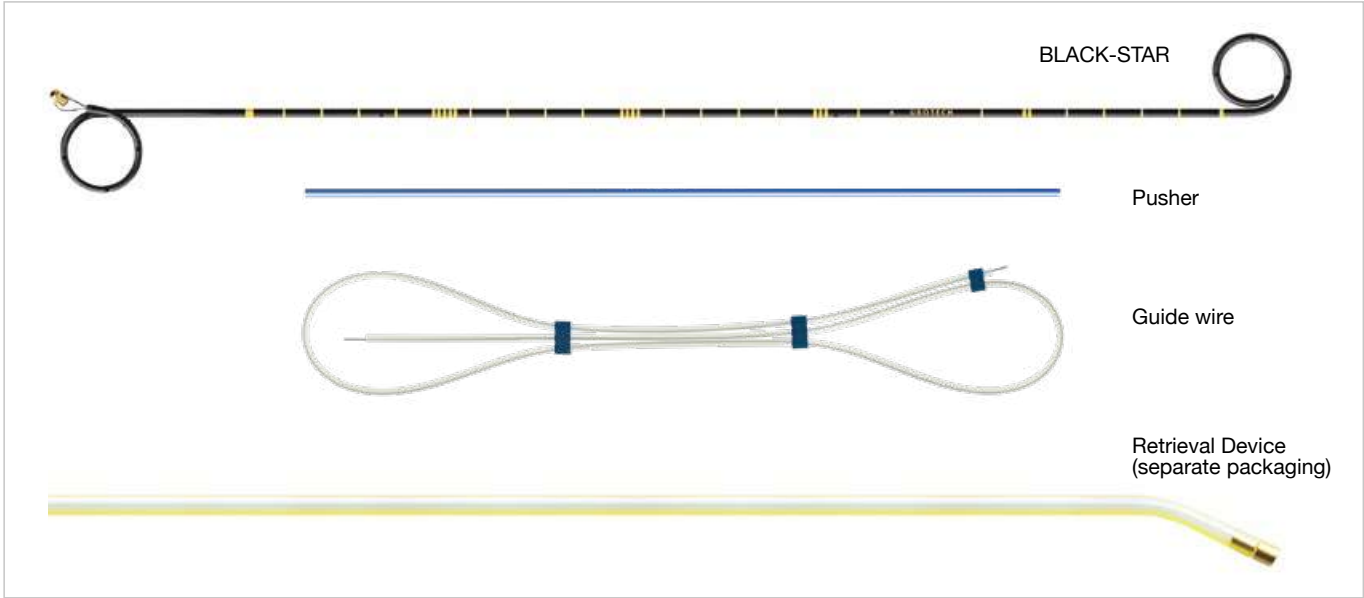

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Drainage openings in pigtails and shaft
- ◆ Markings for exact positioning
- ◆ X-ray visible
- ◆ Bladder sided magnet size S - FR/CH7
- ◆ No need for another Cystoscopy and/or anesthesia to remove the Magnetic Black-Star



Pusher

- ◆ Length: 40 cm

Guide wire standard

- ◆ PTFE-coated

Retrieval Device with magnet

- ◆ Soft PUR
- ◆ Easy tip orientation via marking line
- ◆ Tiemann tip, 30°



URETERAL STENTS

Open tip with Retrieval Device								Packaging unit = 1 Pcs
Article code w/o guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Size Stent-magnet (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
MR-414810	MR-364810	4.8	10	7	145	0.028"	straight	standard
MR-414812	MR-364812	4.8	12	7	145	0.028"	straight	standard
MR-414815	MR-364815	4.8	15	7	145	0.028"	straight	standard
MR-414818	MR-364818	4.8	18	7	145	0.028"	straight	standard
MR-414820	MR-364820	4.8	20	7	145	0.028"	straight	standard
MR-414822	MR-364822	4.8	22	7	145	0.028"	straight	standard

Open tip w/o Retrieval Device								Packaging unit = 1 Pcs
Article code w/o guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Size Stent-magnet (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
MR-404810	MR-304810	4.8	10	7	145	0.028"	straight	standard
MR-404812	MR-304812	4.8	12	7	145	0.028"	straight	standard
MR-404815	MR-304815	4.8	15	7	145	0.028"	straight	standard
MR-404818	MR-304818	4.8	18	7	145	0.028"	straight	standard
MR-404820	MR-304820	4.8	20	7	145	0.028"	straight	standard
MR-404822	MR-304822	4.8	22	7	145	0.028"	straight	standard

Retrieval Device			Packaging unit = 1 Pcs
Article code fits with FR/CH	Size (FR/CH)	Length (cm)	
R-Device-S	9	40	


URETERAL STENTS

MAGNETIC BLACK-STAR® XS and Retrieval Device with magnet

NEW!




SET COMPONENTS

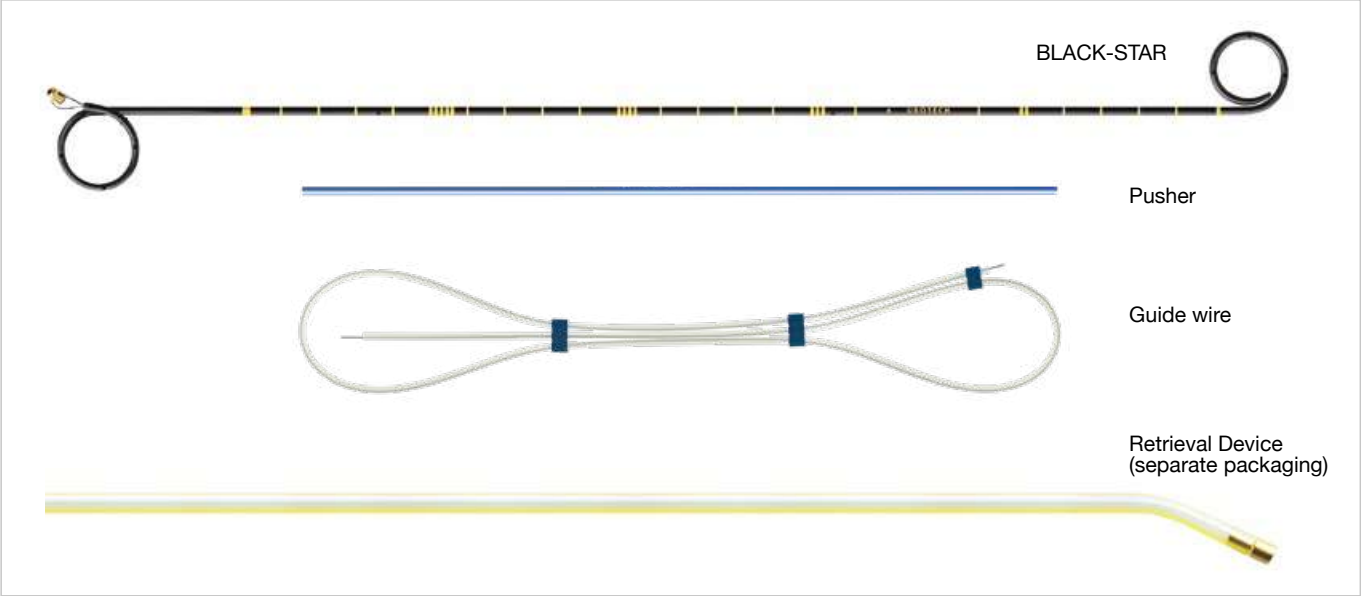

Ureteral Stent

- ◆ Polyurethane (PUR)
- ◆ Drainage openings in pigtails and shaft
- ◆ Markings for exact positioning
- ◆ X-ray visible
- ◆ Bladder sided magnet size XS - FR/CH 4.8
- ◆ Can be used with Cysto-/urethrosopes with a working channel of CH/FR 6 or more
- ◆ No need for another Cystoscopy and/or anesthesia to remove the Magnetic Black-Star XS

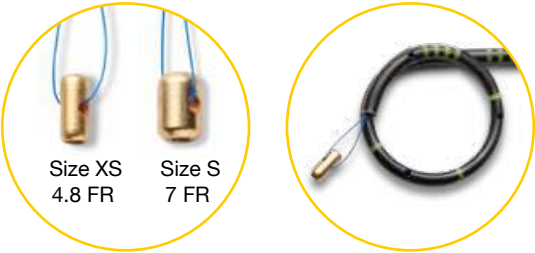

Pusher

- ◆ Length: 40 cm

Guide wire standard

- ◆ PTFE-coated

Retrieval Device with magnet

- ◆ Soft PUR
- ◆ Easy tip orientation via marking line
- ◆ Tiemann tip, 30°



URETERAL STENTS

Open tip with Retrieval Device								Packaging unit = 1 Pcs
Article code w/o Guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Size Stent-magnet (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
MR-414810-XS	MR-364810-XS	4.8	10	4.8	145	0.028"	straight	standard
MR-414812-XS	MR-364812-XS	4.8	12	4.8	145	0.028"	straight	standard
MR-414815-XS	MR-364815-XS	4.8	15	4.8	145	0.028"	straight	standard
MR-414818-XS	MR-364818-XS	4.8	18	4.8	145	0.028"	straight	standard
MR-414820-XS	MR-364820-XS	4.8	20	4.8	145	0.028"	straight	standard
MR-414822-XS	MR-364822-XS	4.8	22	4.8	145	0.028"	straight	standard

Open tip w/o Retrieval Device								Packaging unit = 1 Pcs
Article code w/o Guide wire	Article code with guide wire	Size (FR/CH)	Length (cm)	Size Stent-magnet (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
MR-404810-XS	MR-304810-XS	4.8	10	4.8	145	0.028"	straight	standard
MR-404812-XS	MR-304812-XS	4.8	12	4.8	145	0.028"	straight	standard
MR-404815-XS	MR-304815-XS	4.8	15	4.8	145	0.028"	straight	standard
MR-404818-XS	MR-304818-XS	4.8	18	4.8	145	0.028"	straight	standard
MR-404820-XS	MR-304820-XS	4.8	20	4.8	145	0.028"	straight	standard
MR-404822-XS	MR-304822-XS	4.8	22	4.8	145	0.028"	straight	standard

Retrieval Device			Packaging unit = 1 Pcs
Article code fits with FR/CH	Size (FR/CH)	Length (cm)	
R-Device-S	9	40	

Note:
When inserting the Magnetic Black Star endoscopically, pay attention to the working channel of the endoscope used. This should be at least 1 FR/CH larger than the ureteral stent.


NEPHROSTOMY

RENODRAIN-YELLOW Nephrostomy Puncture- and Exchange Set

Seldinger technique technique, PC-treated, with pigtail




SET COMPONENTS

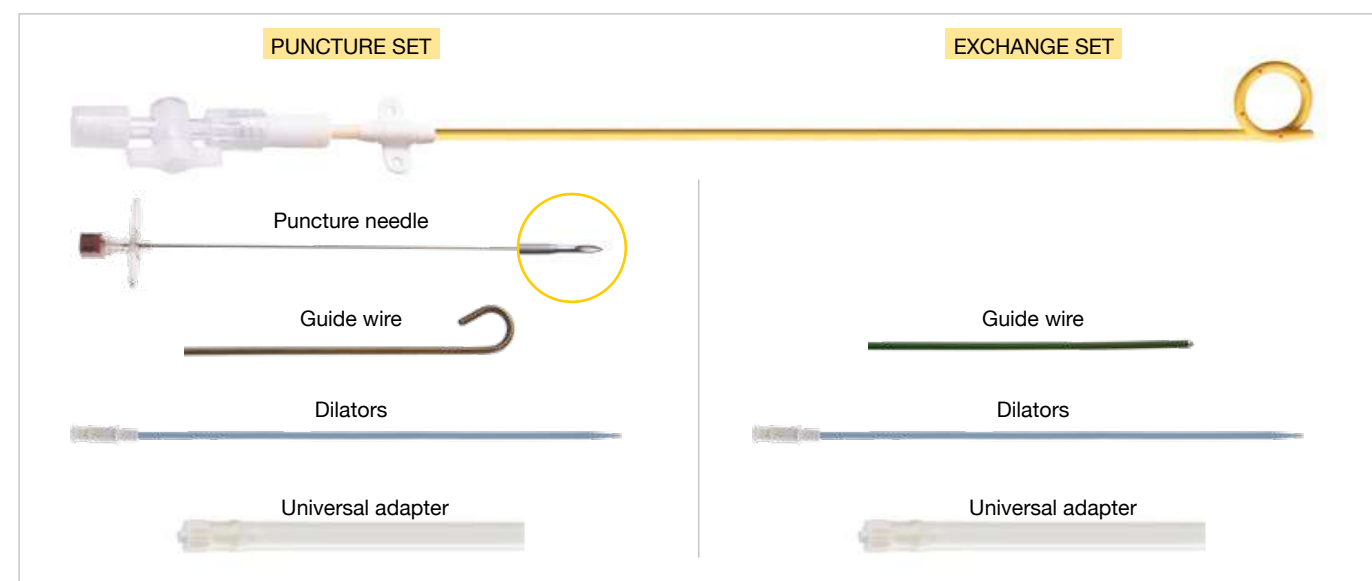

Nephrostomy Catheter

- ◆ Polyurethane (PUR)
- ◆ PC-treated to minimize the risk of infection and encrustation
- ◆ Drainage openings in pigtail
- ◆ Luer-lock Connector

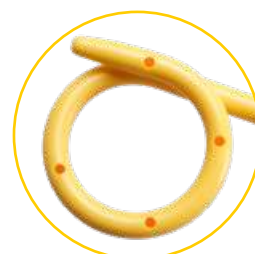

Puncture Needle

- ◆ 2-part needle
- ◆ Medical stainless steel
- ◆ 18 G, length: 20 cm
- ◆ Improved ultrasonic visibility due to sandblasted tip

Guide wire (Lunderquist)

- ◆ PTFE-coated guide wire in puncture set with flexible 1.5mm J-Tip
- ◆ PTFE-coated guide wire in exchange set with flexible straight tip

Universal Adapter

- ◆ Male Adapter with Luer-Lock Connector
- ◆ 10 cm Silicone tube for connecting an urine bag
- ◆ PVC free

NEPHROSTOMY

Single catheter

Packaging unit = 10 Pcs

Article code	Size (FR/CH)	Length (cm)
RP-410640-K	6	24

Puncture Set

Packaging unit = 5 Pcs

Article code	Size (FR/CH)	Length (cm)	Dilators (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip (r=1.5 mm)	Guide wire stiffness
RP-400640-K	6	24	5/6	70	0.025"	J-Tip	hard

Exchange Set

Packaging unit = 5 Pcs


Article code	Size (FR/CH)	Length (cm)	Dilators (FR/CH)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness
RP-400642-K	6	24	5/6	70	0.025"	straight	hard

NEPHROSTOMY

RENODRAIN-INTEGRALE


Nephrostomy Puncture- and Exchange Set, Seldinger technique

SET-KOMPONENTEN

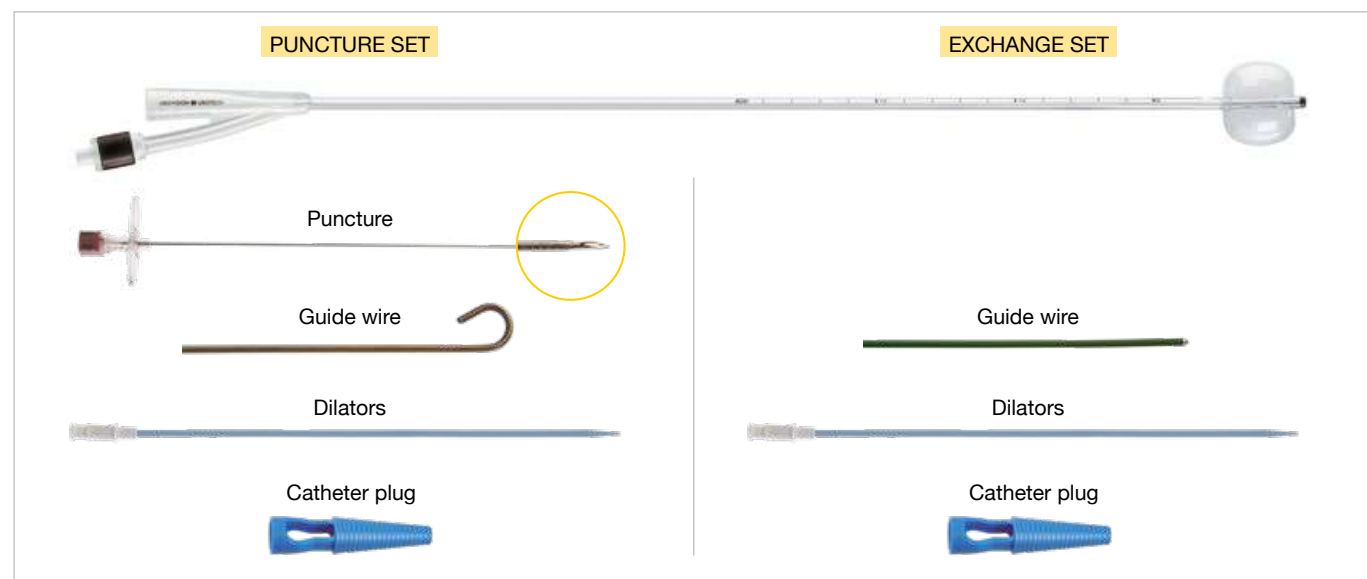

Balloon Catheter

- 100% full silicone catheter with integrated balloon
- Hole behind the balloon for calyx drainage
- X-ray visible contrast stripes
- Markings for exact positioning
- Reinforced funnel attachment for secure connection of the urine bag

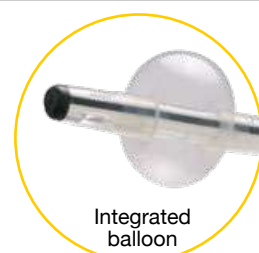

Puncture Needle

- 2-part needle with facet cut
- Medical stainless steel
- 18 G, length: 20 cm
- Improved ultrasonic visibility due to sandblasted tip

Guide wire (Lunderquist)

- PTFE-coated guide wire in puncture set with flexible 1.5mm J-Tip
- PTFE-coated guide wire in exchange set with flexible straight tip

NEPHROSTOMY

Single catheter				Packaging unit = 10 Pcs
Article code	Size (FR/CH)	Length (cm)	Balloon (ml)	
RN-431026	10	30	1.5	
RN-421040	10	40	1.5	

Puncture Set									Packaging unit = 5 Pcs
Article code	Size (FR/CH)	Length (cm)	Balloon (ml)	Dilators (FR/CH) *with Split sheath	Guide wire length (cm)	Guide wire (inch)	Guide wire tip (r=1.5 mm)	Guide wire stiffness	
RN-401040	10	40	1.5	6/8/10/12*	70	0.035"	J-Tip	hard	

Exchange Set									Packaging unit = 5 Pcs
Article code	Size (FR/CH)	Length (cm)	Ballon (ml)	Dilators (FR/CH) *with Split sheath	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness	
RN-411040	10	40	1.5	6/8/10/12*	70	0.035"	straight	hard	

Blocking syringe		PU = 200 Pcs
Article code	(ml)	
SP-2100-10	10	

- Aqua-Dest. glycerin 10%




SUPRAPUBIC BLADDER DRAINAGE

SUPRA-FIX Puncture- and Exchange Set

With pigtail, with funnel

SET COMPONENTS

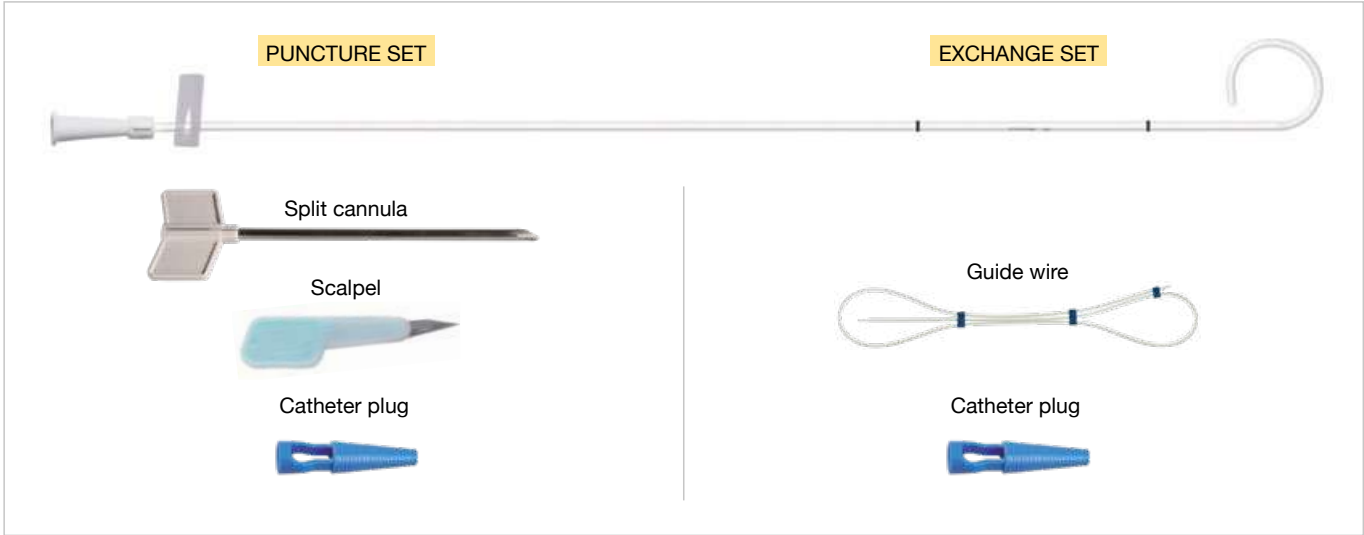

Pigtail catheter

- ◆ Polyurethane (PUR)
- ◆ Flexible pigtail to avoid irritation
- ◆ Funnel for connection with the urine bag
- ◆ Two markings for exact positioning

Split Cannula

- ◆ Medical stainless steel
- ◆ Easy and safe splitting

Guide wire (only in Exchange Set)

- ◆ PTFE-coated

SUPRAPUBIC BLADDER DRAINAGE


Puncture Set					Packaging unit = 5 Pcs
Article code	Size (FR/CH)	Length (cm)	Split cannula (mm)	Split cannula Length (cm)	
SU-081030	8.0	30	3.6	8	
SU-081050	8.0	50	3.6	8	

Exchange Set							Packaging unit = 5 Pcs
Article code	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire (inch)	Guide wire tip	Guide wire stiffness	
SU-081051	8.0	50	100	0.038"	straight	standard	

SUPRAPUBIC BLADDER DRAINAGE

SUPRAVISION INTEGRAL STANDARD
Suprapubic Puncture- and Exchange Set, Seldinger technique

SET COMPONENTS

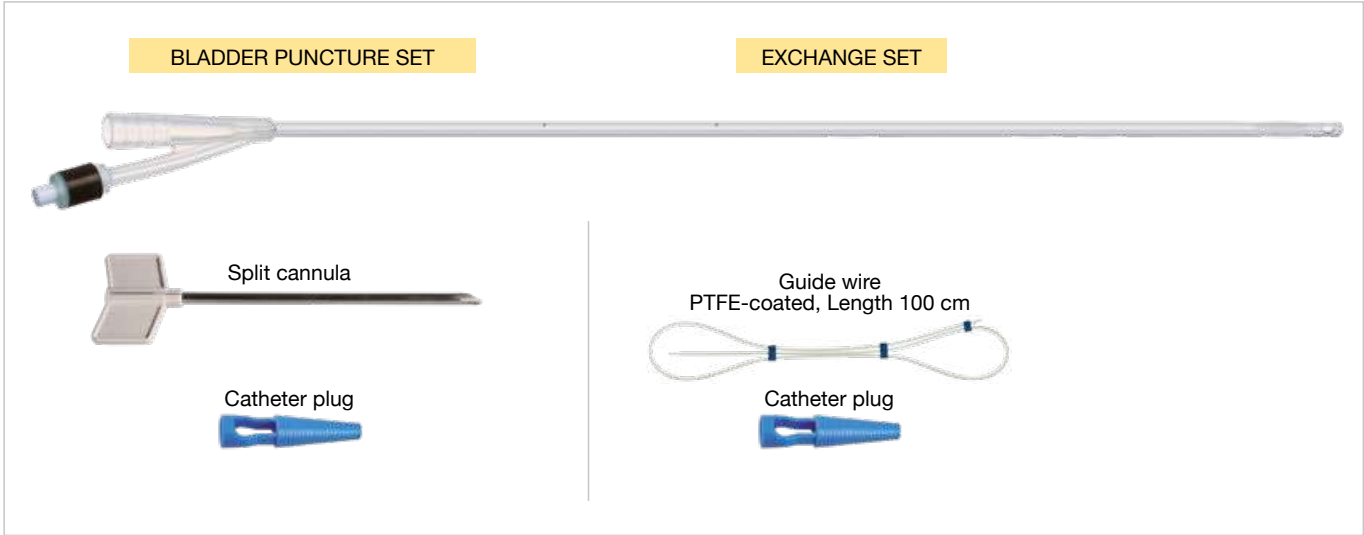

Suprapubic Catheter

- 100% full silicone catheter with integrated balloon
- Central open tip
- Transparent with integrated X-ray contrast strip
- Rounded tip for atraumatic insertion of the catheter
- Reinforced funnel for connection with the urine bag
- Excellent flow-through rate



VISIOSLIDE-Special Coating

- Increased lubricating properties due to special coating
- Wear-resistant inside and outside
- Very low incrustation propensity



SUPRAPUBIC BLADDER DRAINAGE

Puncture Set with VISIOSLIDE special coating						Packaging unit = 5 Pcs
Article code VISIOSLIDE	Size (CH/FR)	Length (cm)	Balloon (ml)	Guide wire (inch)	Split cannula (mm)	Split cannula (cm)
SV-300835-SC	8	35	1.5	0.028"	3.6	8
SV-301040-SC	10	40	3	0.032"	4.6	12

Exchange Set with VISIOSLIDE special coating					Packaging unit = 5 Pcs	
Article code VISIOSLIDE	Size (CH/FR)	Length (cm)	Balloon (ml)	Guide wire (inch)		
SV-330835-SC	8	35	1.5	0.028"		
SV-331040-SC	10	40	3	0.032"		

Single Catheter					Packaging unit = 10 Pcs	
Article code VISIOSLIDE	Size (CH/FR)	Length (cm)	Balloon (ml)	Guide wire (inch)		
SV-310835-SC	8	35	1.5	0.028"		
SV-311040-SC	10	40	3	0.032"		

Blocking syringe		PU = 200 Pcs
Article code	(ml)	
SP-2100-10	10	

- Aqua-Dest. glycerin 10%




TRANSURETHRAL BLADDER DRAINAGE

2-WAY BALLOON CATHETER

Standard and Supersoft

ESCRPTION

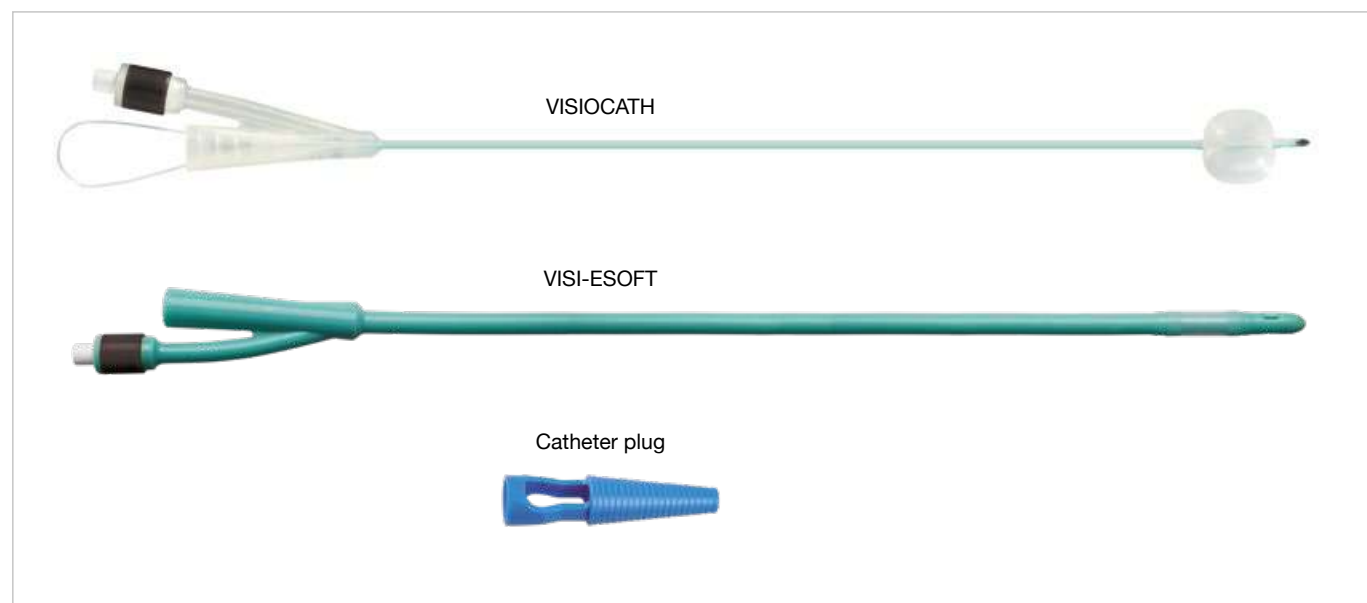

VISIOCATH 2-way balloon catheter

- 100% full silicone for optimal patient comfort
- Good drainage performance due to large inner diameter
- Very soft surface for pleasant feeling on the skin
- With internal plastic stylet
- Very good lubricating properties
- Wear-resistant inside and outside
- Low incrustation propensity

With
special
VISIOSLIDE
coating

VISI-ESOFT 2-way balloon catheter

- 100% full silicone for optimal patient comfort
- Very soft balloon catheter for a very high comfort level for the patient
- Good drainage performance due to large inner diameter
- Reinforced funnel for connection with the urine bag

TRANSURETHRAL BLADDER DRAINAGE

VISIOCATH 2-way balloon catheter

Packaging unit = 10 Pcs

Article code	Size (CH/FR)	Length (cm)	Balloon (ml)
KW-450603-SC	6	30	1.5
KW-450803-SC	8	30	5
KW-451003-SC	10	30	5

VISI-ESOFT 2-way balloon catheter

Packaging unit = 10 Pcs

Article code	Size (CH/FR)	Length (cm)	Balloon (ml)
KG-300810	8	30	3
KG-301010	10	30	5

Blocking syringe

PU = 200 Pcs

Article code	(ml)
SP-2100-10	10

- Aqua-Dest.-glycerin 10 %



SPECIAL PRODUCTS

PYELOPLASTY CATHETERS

For renal pelvis plastic

SET COMPONENTS



Pyeloplasty Catheter

- ◆ Polyurethane (PUR)
- ◆ Drainage openings in pigtail and ureter section
- ◆ Markings for exact positioning
- ◆ X-ray visible

Needle adapter with funnel

- ◆ For connection with the urine bag

Article code	Size (FR/CH)	Ureter-/external part (cm)	Max Guide wire	Pigtail Ø (mm)
EK-620528	4.7	6 / 20	0.028"	15
EK-620628	6	6 / 20	0.035"	15
UR-910435	4	7 / 21	0.025"	15
UR-910535	5	7 / 21	0.028"	15


Packaging unit = 10 Pcs

SPECIAL PRODUCTS

LINZER CATHETER SET

Hypospadias and Epispadias

SET COMPONENTS

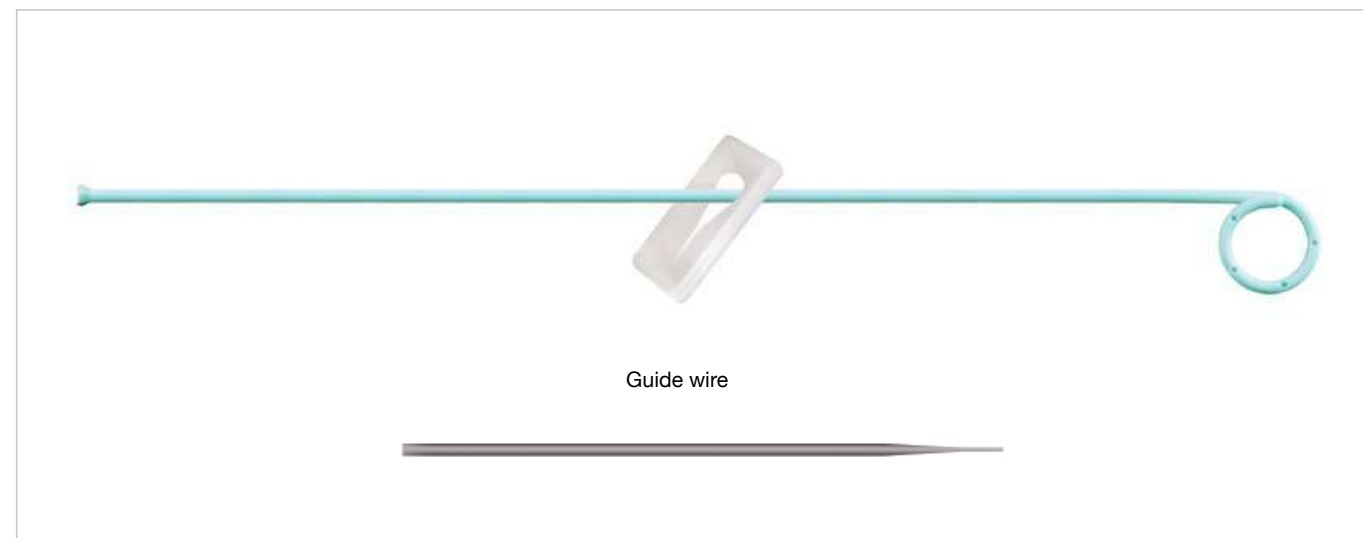

Linzer Catheter set

- ◆ Soft polyurethane (PUR) for best patient comfort
- ◆ Splinting of the urethra after surgical interventions for hypospadias and epispadias
- ◆ Ensuring postoperative bladder drainage
- ◆ Stabilizes the urethra, without suture

Guide wire

- ◆ PTFE-coated

Closed tip						Packaging unit = 10 Pcs
Article code	Size (FR/CH)	Length (cm)	Guide wire length (cm)	Guide wire	Guide wire tip	Guide wire stiffness
EG-350822	8	22	40	0.038"	straight	hard

ORDERING INFORMATION

Order via Fax



+49 8032 973-211

Order via phone



+49 8032 973-210

Order via Email



info@urotech.com

UROVISION GmbH

UROTECH GmbH

Medi-Globe-Str. 1-5

D-83101 Rohrdorf OT Achenmühle

Telefon: +49 (0) 8032 - 973-210

Fax: +49 (0) 8032 - 973-211

E-Mail: info@urotech.com

www.urovision-urotech.de



UROVISION•UROTECH

A MEDI-GLOBE COMPANY

LL-0421 Vers. 1.0 03/2021